

Indications

- Wrinkles and fine lines
- Skin discoloration
- Sun damage skin
- Skin tone and texture
- Tightening of the upper eyelids
- Age spots
- Stretch marks / Striae

Treatments

- Spacing treatments sessions - 1 month apart
- Treatment program : 3 - 4 sessions
- After treatments skin slightly swollen and red for 3 to 4 hours
- Patients experience mild erythema and heat sensation for 2 to 3 hours
- Downtime from 0 to 2 days depending on treatment intensity
- Micro crusting appears the next day, disappears after about 3 days
 - Improvement in skin texture, younger look, wrinkle reduction
 - Better definition of eyelids
 - High patient satisfaction, all skin type

ADVANCED TECHNOLOGY FOR FRACTIONAL SKIN REJUVENATION



Acne

- Female, 22
- Follow-up 30 days 4 Tixel Tx
- No Downtime
- Without topical anesthesia



Facial Treatment

- Female, 66
- Follow-up 7 months 4 Tixel Tx
- 1 day downtime
- Without topical anesthesia



Technical Specifications

Energy Type:	Direct heating
Titanium Tip:	Yes
Tip Coverage Area:	1 cm ² (Facial Tip), 0.3 cm ² (periorbital Tip)
System Weight (w/o HP):	6 kg
Handpieces Weight:	330 grams
Console Dimensions:	L 26 x W 44 x H 30 (cm)
Handpieces:	2 parallel
Settings:	5 -18 ms, 100 - 1000 μm

Better For Patient

- Low pain
- No charring
- Short downtime
- Treat all skin types


Better For Clinic

- Excellent results
- Safe
- Simple and fast
- Small

Exclusive Distributor:

REDO
MARKETING INDONESIA
Your Aesthetic Business Solutions Provider

Jalur Sutera Timur, Ruko Renata Blok 1A/9
Alam Sutera, Tangerang 15144
Tel: 021-2966 1469 Fax: 021-3043 9750
Website: www.redo.co.id Email: enquiry@redo.co.id
All Right Reserved by PT. Redo Marketing Indonesia

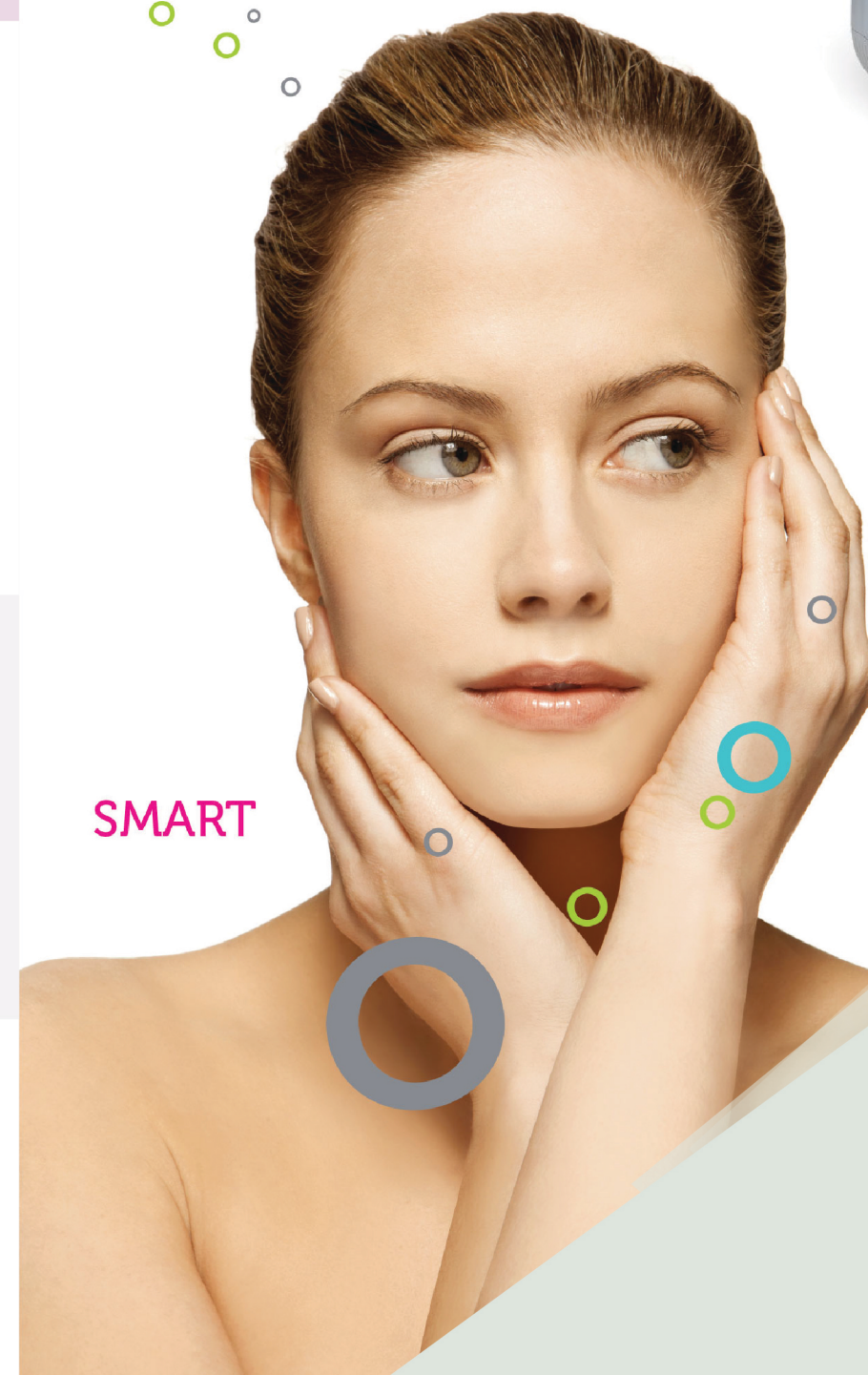


tixel[®]
2



What Is Tixel® ?

- A novel non-ablative & ablative skin rejuvenation system powered by Novoxel's patented titanium tip technology.
- Tixel treats delicate facial skin, including periorbital, eyelids, the neck and décolleté safely and quickly with low pain and 0-2 days downtime.
- Tixel employs direct heat and is radiation-free.

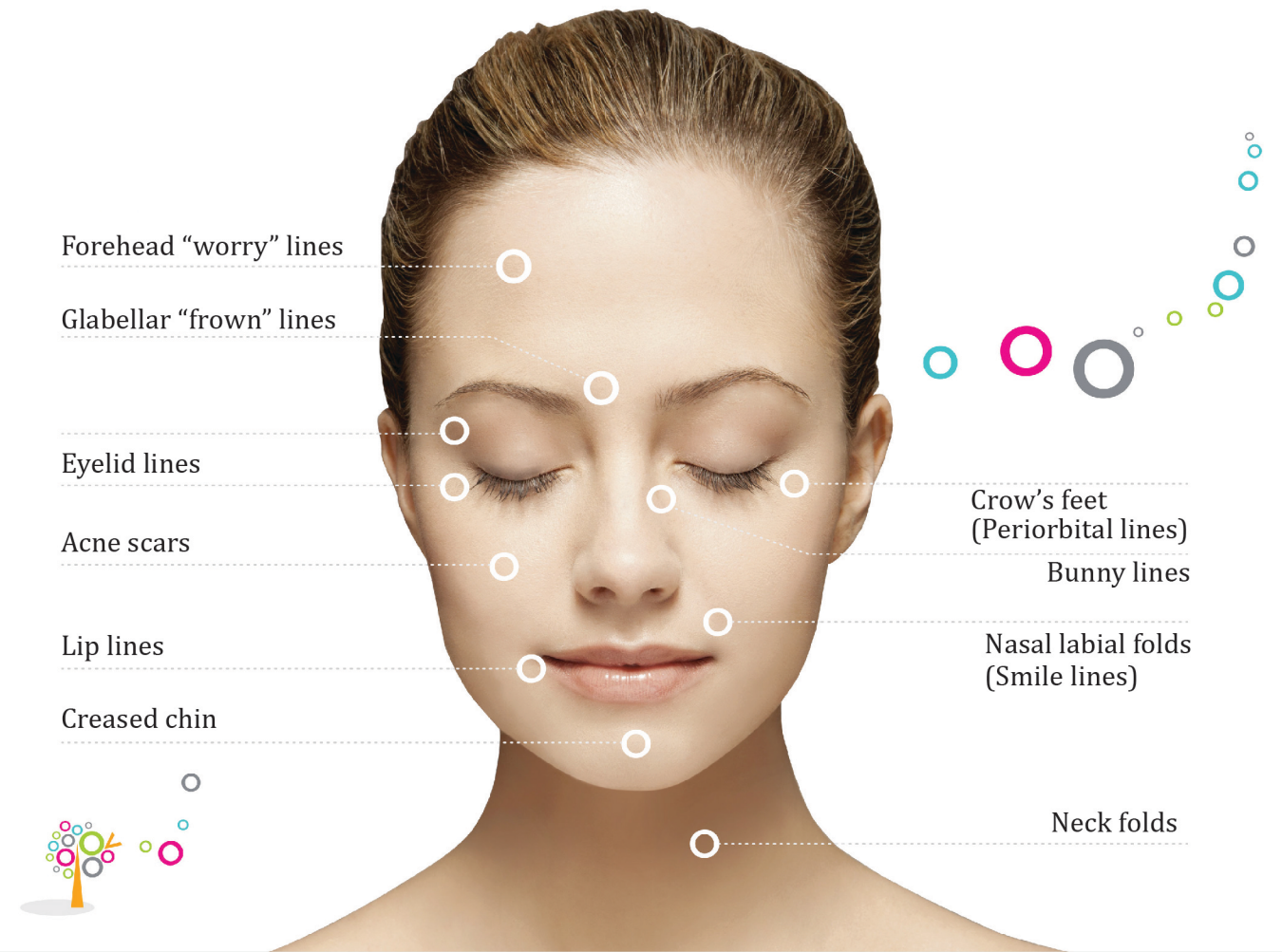

Our Benefits



MULTIPLE TREATMENTS

Skin rejuvenation, crow's feet, open pores, acne, stretch marks....



MULTIPLE AREAS

Face, neck, décolleté hands, arms & body



LOW PAIN

Analgesic creams are not required



STERILE TITANIUM TIP

The Tip is biocompatible and sterile during treatment



SMALL & PORTABLE


Small, light, easy to carry



RADIATION FREE

Tixel employes direct heat

Patented - Scientifically Validated - Confirmed By Results

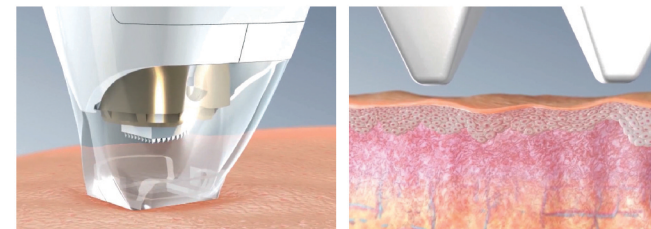

Titanium Tip

- Durable
- Clean
- Sterile
- Precise

Thermal energy is transferred to the skin by a hot titanium Tip, which contacts the skin briefly. There is no bleeding.

The extent of energy transfer is determined by the Tip's exposure time with the skin. Heat transfer is by contact. The Tip does not perforate the tissue.



Periorbital Handpiece



Facial Handpiece